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Instruction Manual *Please read carefully!*

EXPLANATION OF SYMBOLS

The following signs and symbols on items, packaging and in the instruction manual bear important information:

 This instruction manual belongs to this device. It contains important information about starting up and operation. Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.

 **WARNING**
These warning notes must be observed to prevent any injury to the user.

 **CAUTION**
These notes must be observed to prevent any damage to the device.







 **NOTE**
These notes give you useful additional information on the installation or operation.

 These lancets for blood glucose tests correspond to the requirements of the EU guideline "93/42/EEC of the Council Directive dated 14 June 1993 concerning medical devices" and bear the CE marking (conformance label) "CE 0297".

	Do not re-use		LOT number
	Use by		Catalogue number
	Storage temperature limitation		Manufacturer
			Sterilization by radiation

  **WARNING**
Never use a lancet or a lancing device that belongs to another person. Use a new sterile lancet for each test. Lancets are intended to be used only once. Used lancets are considered hazardous, biological non-biodegradable waste. Dispose of them taking into account that they are capable of spreading infection. Dispose of the lancets so that there is no risk of injury or infection to other people.

USE

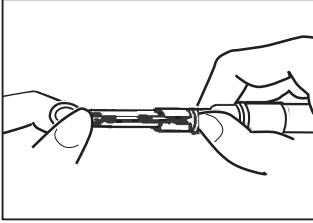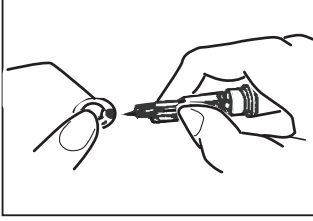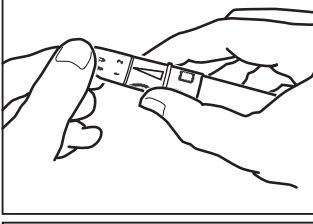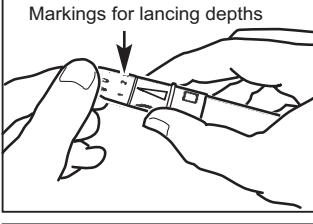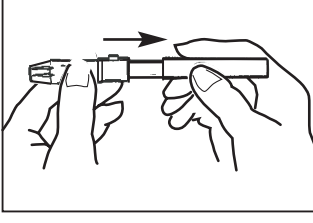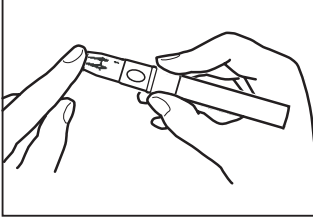- 
 - Wash your hands with soap and warm water. Rinse and dry thoroughly. If necessary, wipe the skin where the blood sample will be taken with an alcohol pad.
- 
 - Open the lancing device by twisting the cap in a clockwise direction and then remove it. Insert the lancet all the way (without turning it) into the cup in the lancing device.
- 
 - Twist off the protective cap of the lancet carefully. Retain this cap to aid disposal of the lancet after use.
- 
 - Replace the protective cap on the lancing device and tighten by turning it in an anti-clockwise direction. To use the AST cap, consult the instruction manual.
- 
 - The lancing device can be set to the sensitivity of your skin. You can adjust the tip to 5 different lancing depths. Twist the cap end in the appropriate direction until the arrow is pointing to the number for the lancing depth you want. Consult the following for the suitable lancing depth:
 - **1 - 2** for soft or thin skin
 - **3** for normal skin
 - **4 - 5** for thick or callous skin
- 
 - Set the lancing device by extending it until it clicks into place. If it does not click, it is probably already in position from when you inserted the lancet. The lancing device is ready to use. Do not lance your finger before the monitor and test strips are ready to use.
- 
 - Once the monitor and test strips are ready and the area of skin has been cleaned, place the lancing device on a finger tip (preferably at the side) and press the trigger button. Draw a drop of blood by gently massaging the area. Make sure that the blood droplet does not smudge and perform the test as described in the instructions for use.

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MEDISANA AG